

AVAIL

Miracle-Ear**AVAIL**[™]
PRODUCT PORTFOLIO

NEW!

Miracle-Ear**AVAIL**™



WHAT'S INSIDE...

CONTENTS:

Discover Miracle-Ear**AVAIL**™, our NEW line of innovative hearing aids.



Products available in solution levels 5 to 2*



Versatile accessories for a seamless connection

Miracle-Ear AVAIL ™ TECHNOLOGY	4
Miracle-Ear ENERGY ™ ITC R EV	7
Miracle-Ear ENERGY ™ RIC R EV	8
Miracle-Ear ENERGY ™ BTE R EV	9
Miracle-Ear MINI ™ IIC 10 EV	11
Miracle-Ear MINI ™ CIC 10 EV	12
Miracle-Ear MINI ™ CIC 312 EV	13
Miracle-Ear CONNECT ™ BTE SP 13 EV	15
Miracle-Ear EASY ™ BTE 13 SP EV	17
Miracle-Ear® App	18
Miracle-Ear AVAIL ™ Accessories Overview	19
LockFit™ RIC fitting solutions	20
LockFit™ BTE fitting solutions	23

*Availability varies by product type

Miracle-EarAVAIL™ TECHNOLOGY

Seamless innovation.

For more than 70 years, Miracle-Ear® has helped millions of people live better, more fulfilling lives through personalized expert care and by providing innovative hearing aid solutions.

Miracle-EarAVAIL™ elevates our product portfolio with our most advanced approach to natural sound quality, custom craftsmanship, seamless personalization, and enhanced connection to effortlessly fit into our customer's lives.

NEW AUDIOLOGICAL FEATURES



M-E Mode™¹

Provides instant optimized adjustments and puts the power of artificial intelligence at the customer's fingertips for next-generation sound performance in the most challenging listening environments.



Custom Tinnitus +™¹

Tool allows the hearing care professional to set personalized shaping of the tinnitus stimulus for the customer.



24 Band Technology¹

Expanded options to personalize for harder to fit customers.



Direct Streaming

Available for iOS® and Android™ in every product family.

M-E Mode™ Key Features:



Real-time adjustments made for the listener



Reduces listening efforts



Same or better than optimized hearing aid memory



Setup "double-tap" when fitting the hearing aid for instant adjustments

SEAMLESS INTELLIGENCE.

Our latest technology provides superior sound quality with advanced features and capabilities.

Miracle-EarAVAIL™				
	Level 5	Level 4	Level 3	Level 2
M-E Mode™ with Tap Control	•			
Custom Tinnitus +™	•			
24 Band Personalization	•			
Hands-Free with 2-way Audio (iOS®)	•	•	•	•
Auto On/Off	•	•	•	•
Tap Control User Settings	•	•	•	•
REM Perfect™ Fitting	•	•	•	•

¹Available in Level 5 only.

Miracle-EarAVAIL™ PERSONALIZATION

Seamless craftsmanship.




Miracle-EarAVAIL™ hearing aids are expertly designed in the United States, handcrafted and personalized for your customer's unique ear to provide exceptional sound quality.

Personalization and custom craftsmanship have led to a seamless hearing experience. Customers now have a range of hearing aid style options that fit into their unique lifestyles, including our custom ITE hearing aids and molds handcrafted for the most comfortable fit.

From our smallest to our most powerful, customers will experience realistic and genuine sound quality in every listening environment without the need to do anything extra.



Miracle-Ear**ENERGY**,™ rechargeable hearing aids

Solution Level	ITC	RIC	BTE
5	MEENERGY™ ITC R EV	MEENERGY™ RIC R EV	MEENERGY™ BTE R EV
4			
3	Miracle-Ear AVAIL ™ Technology	Miracle-Ear AVAIL ™ Technology	Miracle-Ear AVAIL ™ Technology
2			

NEW!
MEENERGY™
ITC R EV



KEY FEATURES

- Small and recessed for the utmost discretion
- Artificial Intelligence and Embedded Sensor Technology
- Push button user control
- Stream from Apple iOS® and select Android™ devices
- Wireless connectivity
- 23 Hours of superior hearing, including 4 hours of streaming
- Tinnitus Technology
- Compatible with all MELINK™ Accessories




Compatible with Miracle-Ear® App only



Integrated sensor for double tap control



Rechargeable

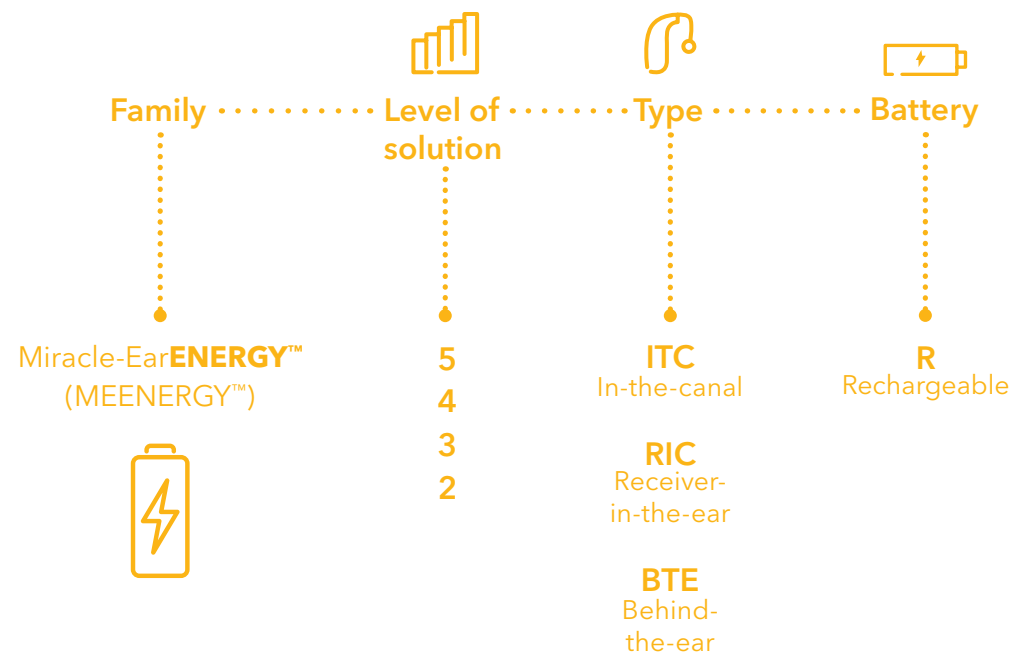


Solution levels

5
4
3



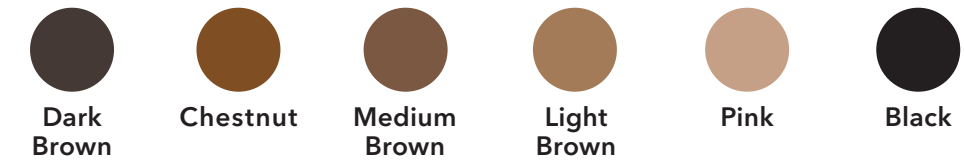

Color Guide

Faceplate Colors



Shell Colors



Model	MEENERGY™ ITC R EV			
Receiver type				
Performance	Miracle-Ear AVAIL ™	Miracle-Ear AVAIL ™	Miracle-Ear AVAIL ™	Miracle-Ear AVAIL ™
Matrix (2 ccm)	110/40 dB	115/50 dB	120/60 dB	130/70 dB
Battery	Lithium-ion Encased 312			
Size	Lithium-ion Encased 312	Lithium-ion Encased 312	Lithium-ion Encased 312	Lithium-ion Encased 312

NEW!
MEENERGY™
RIC R EV



Shown in Bright White with Sterling

KEY FEATURES

- Artificial Intelligence and Embedded Sensor Technology
- Stream from Apple iOS® and select Android™ devices
- Wireless connectivity
- 23 Hours of superior hearing, including 4 hours of streaming
- Tinnitus Technology
- Telecoil
- Compatible with MECROS EV
- Compatible with all MELINK™ accessories



Solution levels

- 5
- 4
- 3
- 2

Color Guide

Standard Colors




NEW!
MEENERGY™
BTE R EV



Shown in Black

KEY FEATURES

- Artificial Intelligence and Embedded Sensor Technology
- Stream from Apple iOS® and select Android™ devices
- Wireless connectivity
- 23 Hours of superior hearing, including 4 hours of streaming
- Tinnitus Technology
- Telecoil
- Compatible with all MELINK™ accessories

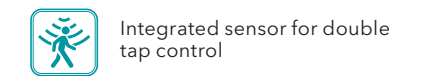

Solution levels

- 5
- 4
- 3
- 2

Color Guide

Standard Colors



Model	MEENERGY™ RIC R EV			
Type	Miracle-Ear ^{AVAIL} ™			
Performance	Matrix (2 ccm)			
	107/40 dB	115/50 dB	120/60 dB	130/70 dB*
Battery	Lithium-ion Encased 312			
Size	Lithium-ion Encased 312	Lithium-ion Encased 312	Lithium-ion Encased 312	Lithium-ion Encased 312

*Requires AP receiver.

Model	MEENERGY™ BTE R EV
Type	Miracle-Ear ^{AVAIL} ™
Performance	Matrix (2 ccm)
	130/70 dB
Battery	Lithium-ion Encased 312
Size	Lithium-ion Encased 312

Miracle-Ear**MINI**,™ our smallest hearing aids


Solution Level	IIC 10	CIC 10	CIC 312
5	MEMINI™ IIC 10 EV	MEMINI™ CIC 10 EV	MEMINI™ CIC 312 EV
4			
3	Miracle-Ear AVAIL ™ Technology	Miracle-Ear AVAIL ™ Technology	Miracle-Ear AVAIL ™ Technology
2			

NEW!
MEMINI™
IIC 10 EV



Shown in Black

- KEY FEATURES**
- Natural sound quality and speech understanding in noise
 - For individual handling and fitting needs
 - Tinnitus Technology
 - 4-7 day battery life

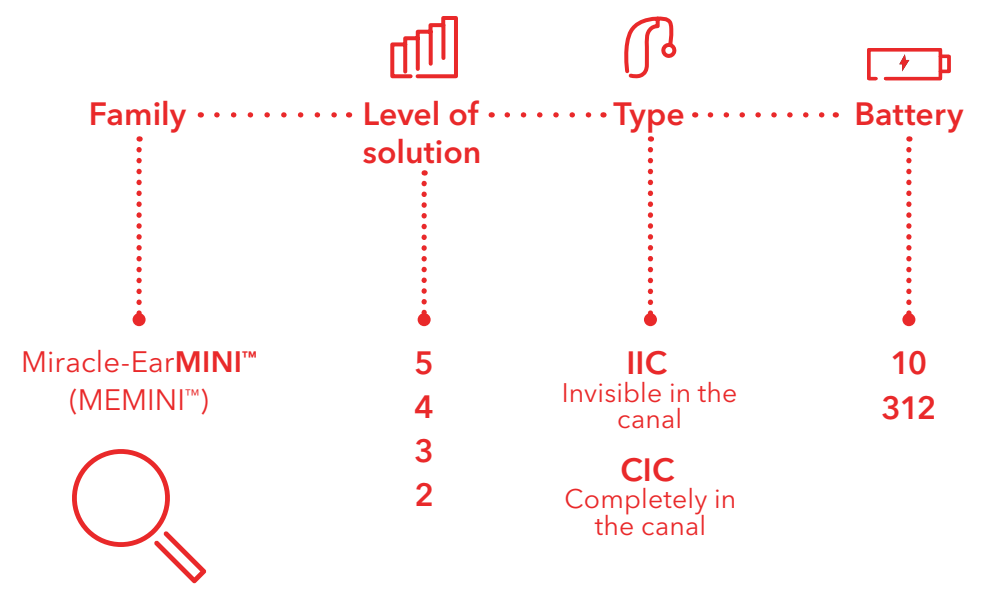




Color Guide

Solution levels

5	Dark Brown	Medium Brown	Light Brown	Pink	Black
4					
3					

Shell Colors

Black

Model	MEMINI™ IIC 10 EV		
Type			
Performance	Miracle-Ear AVAIL ™	Miracle-Ear AVAIL ™	Miracle-Ear AVAIL ™
Matrix (2 ccm)	110/35 dB	110/40 dB	115/50 dB
Battery			
Size	10	10	10

NEW!
MEMINI™
CIC 10 EV



Shown in Light Brown and Clear

KEY FEATURES

- Natural sound quality and speech understanding in noise
- Tinnitus Technology
- Optional push button user control
- 4-7 day battery life

Solution levels

- 5
- 4
- 3

Color Guide


Faceplate Colors



Shell Colors



NEW!
MEMINI™
CIC 312 EV

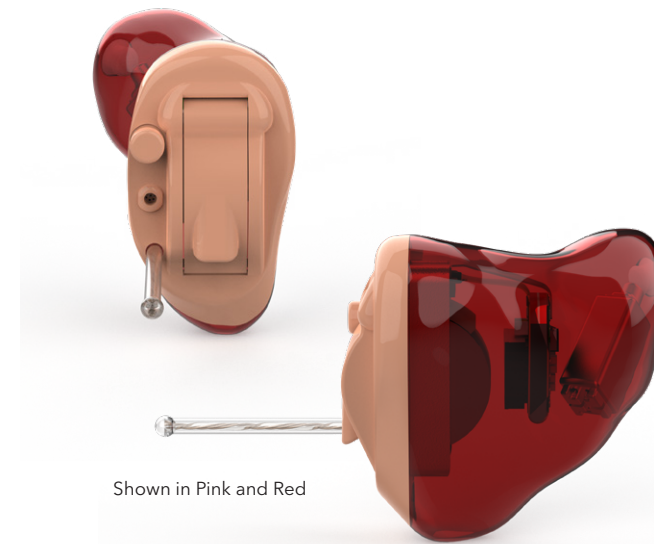

Shown in Pink and Red

KEY FEATURES

- Natural sound quality and speech understanding in noise
- Artificial Intelligence and Embedded Sensor Technology
- Revolutionary antenna design built into the retrieval line
- Stream from Apple iOS® and select Android™ devices
- Wireless connectivity
- Tinnitus Technology
- Optional push button user control
- Compatible with all MELINK™ Accessories
- 4-7 day battery life



Compatible with Miracle-Ear® App only



Integrated sensor for double tap control

Solution levels

- 5
- 4
- 3
- 2

Color Guide

Faceplate Colors



Shell Colors



Model	MEMINI™ CIC 10 EV			
Type				
Performance	Miracle-EarAVAIL™			
Matrix (2 ccm)	110/40 dB	115/50 dB	120/60 dB	130/70 dB
Battery	Size			
	10	10	10	10

Model	MEMINI™ CIC 312 EV			
Type				
Performance	Miracle-EarAVAIL™			
Matrix (2 ccm)	110/40 dB	115/50 dB	120/60 dB	130/70 dB
Battery	Size			
	312	312	312	312

Miracle-EarCONNECT™ Bluetooth® hearing aids




Solution Level	BTE 13
5	MECONNECT™ BTE SP 13 EV
4	
3	Miracle-EarAVAIL™ Technology


NEW!
MECONNECT™
BTE SP 13 EV



KEY FEATURES

- Artificial intelligence and embedded sensor technology
- Tinnitus Technology
- Stream from Apple iOS® and select Android™ devices
- Wireless connectivity
- Artificial intelligence and embedded sensor technology
- Telecoil
- 6-9 day battery life
- Compatible with all MELINK™ Accessories


-  Compatible with Miracle-Ear® App only
-  Integrated sensor for double tap control
-  Advanced acoustic sensors

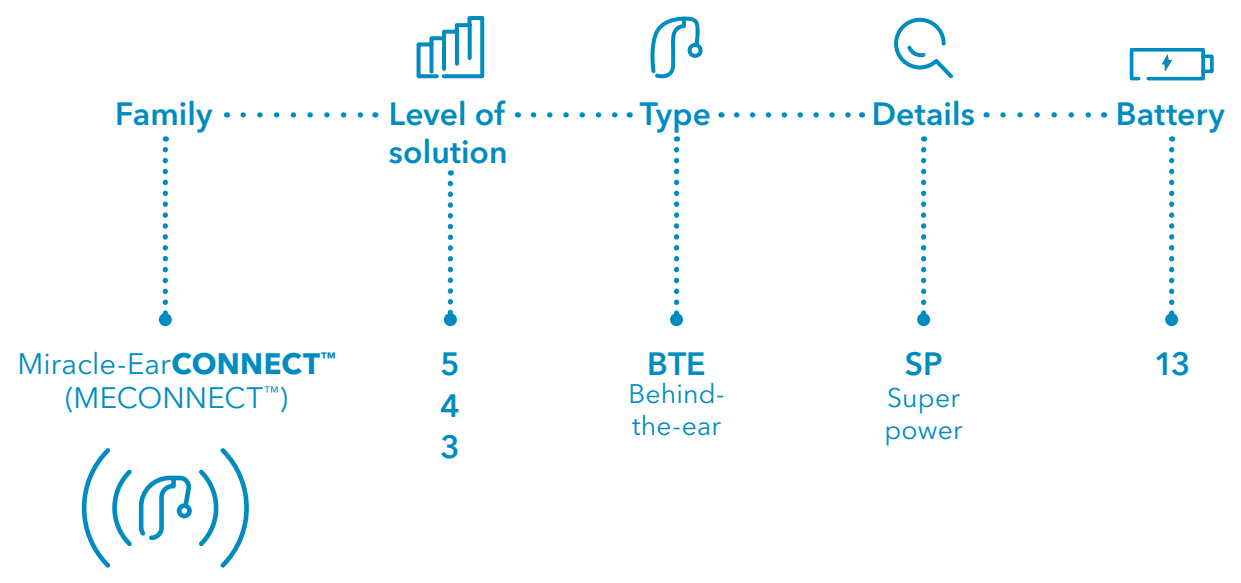


- Solution levels
- 5
 - 4
 - 3

Color Guide

Standard Colors

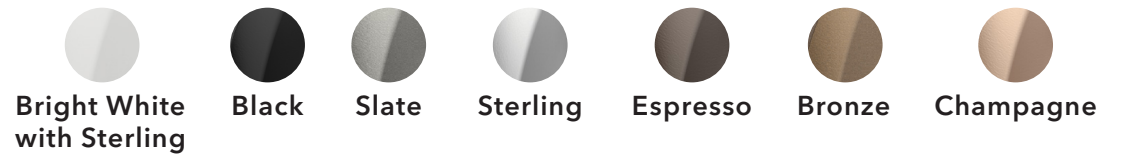

Model	MECONNECT™ BTE SP 13 EV
Receiver type	
Performance	Miracle-EarAVAIL™
Matrix (2 ccm)	140/80 dB
Battery	
Size	13

Miracle-Ear**EASY**,™ everyday features

Solution Level

BTE 13

MEEASY™ BTE SP 13 EV



Miracle-Ear**AVAIL**™ Technology

2



NEW!
MEEASY™
BTE SP 13 EV




Shown in Sterling

KEY FEATURES

- Artificial Intelligence and Embedded Sensor Technology
- Stream from Apple iOS® and select Android™ devices
- Wireless connectivity
- Tinnitus Technology
- Telecoil
- Compatible with all MELINK™ Accessories
- 6-9 day battery life

-  Compatible with Miracle-Ear® App only
-  Integrated sensor for double tap control

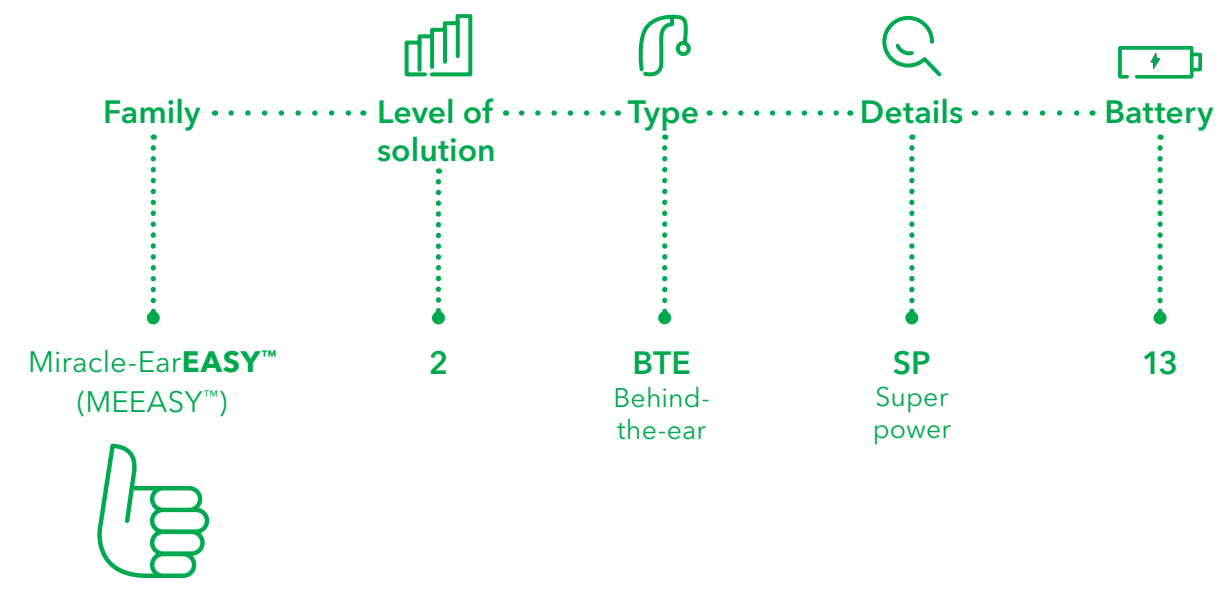



Solution levels

2

Color Guide

Standard Colors

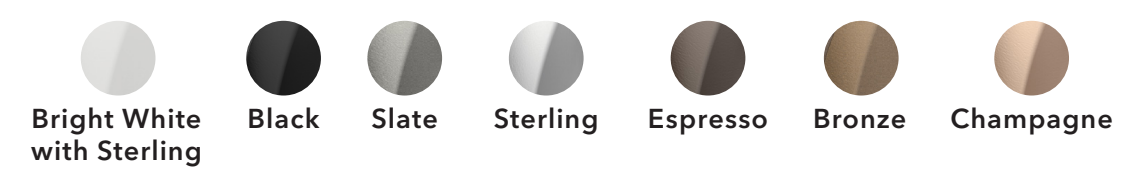

Model	MEEASY™ BTE SP 13 EV
Receiver type	 Miracle-Ear AVAIL ™
Performance	Matrix (2 ccm) 140/80 dB
Battery	Size 13

Miracle-Ear® App


Give your customers full control of their hearing by fusing superb sound with enhanced connectivity.

Our revolutionary app launches Miracle-Ear® into the future. The Miracle-Ear App allows wearers to book an appointment, search locations via the store locator and access the knowledge center to get the latest articles and news, guidance on hearing aids, tutorials and more! Once paired with Miracle-Ear^{AVAIL}™ Bluetooth® hearing aids, wearers can adjust hearing aid volume and swipe to change the program to optimize their hearing for almost every situation.

Visit the Miracle-Ear App Support website accessible in the Operations Manual via the Partner Portal to learn more about the App. The support website allows you to:

- Tour the app to understand features and functionality
- Check compatibility between smartphone and hearing aid models
- Review frequently asked questions
- Review the Miracle-Ear App User Manual

If you or your customer need additional assistance, our Support Team is available to help. Submit a ticket through Partner Portal or give them a call (888) 386-3607.

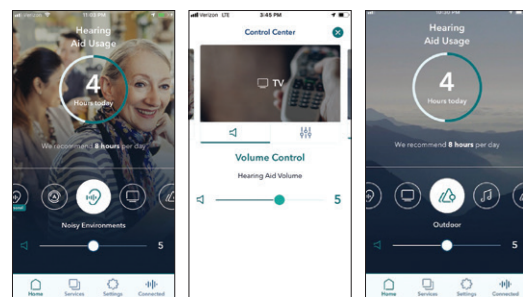

Download for free from the Apple App Store or GooglePlay Store.



**See legal note on final page in this document.

Miracle-Ear^{AVAIL}™ Accessories Overview

Rechargeability



MECHARGE™ start RIC/BTE AV & custom

Offers standard table-top charging that lets your customer know when hearing aids are charged.



MECHARGE™ fullpower RIC/BTE AV & custom

Offers power bank charging on-the-go with up to three charges. This charger lets customers know when their rechargeable hearing aids are charged, automatically turning the devices on when removed from the charger.



MECHARGE™ fullpower RIC/BTE T AV

In just seven minutes, our pocket-sized charger provides 35 hours of power and can completely charge RIC R or BTE R hearing aids up to four times.

Remote Controls & Microphones



MELINK™ Remote AV

Easily control memories and volume, mute hearing aids and turn other special features on and off.



MELINK™ Remote Mic AV

Customers can enjoy one-on-one conversations in noisy environments by attaching it to a conversation partner's clothing.



MELINK™ Remote Mic+ AV

Enables easy audio streaming and one-on-one conversations in noisy environments.

Audio Streamers



MELINK™ TV Streamer AV

Stream audio from your TV or other electronic audio source directly to your hearing aids. It offers excellent sound quality, is easy to use and supports both analog and digital input sources.

NEW!
LockFit™

Compatible hearing aids

MEENERGY™ RIC R EV

RIC fitting solutions

RIC earmolds provide a distinct advantage in that the receiver can be changed out in the office. In addition, our embedded AP receiver allows us to disassemble the receiver configuration and re-orient the components, fitting the customer comfortably.

RIC Earmolds / Digital SLS / HS Silicone 60 / HS Silicone 40

FIND ME: Shop > Earmolds > RIC Earmold

Tubing Lengths:


1, 2, 3, 4, 5

Gain:


40 dB, 50 dB, 60 dB

Receivers:

SnapFit
(six-pin connector)

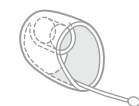

RIC
Hollow Canal

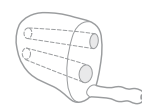

RIC
Canal Solid



RIC
Canal Lock



RIC
Skeleton



RIC
Helix



RIC
Helix Occluded

Absolute Power (AP) Receivers / Digital SLS

Only available in clear

FIND ME: Shop > Absolute Power (AP) Receivers

See specific models on prosite

Tubing Lengths:

1, 2, 3, 4, 5

Gain:

40 dB, 50 dB, 60 dB, 70 dB

Receivers:

SnapFit
(six-pin connector)



RIC AP
Canal



RIC AP
Canal Lock

Pro tip: When measuring the receiver wire for an AP receiver, report the number that correlates with the top of the entrance to the ear canal.

Stock solutions



Open
Comfort Bud
5 mm, 6 mm, 7 mm,
9 mm, 11 mm



Occluded
Comfort Bud
4 mm, 5 mm, 6 mm,
7 mm, 9 mm, 11 mm



Power Domes
8-10 mm, 10-12 mm,
12-14 mm



Stock Helix
RIC Mold



Stock Hollow
Silicone RIC Mold



Stock
RIC Mold





NEW!
LockFit™

Compatible hearing aids

MEENERGY™ BTE R EV
MECONNECT™ BTE SP 13 EV
MEEASY™ BTE SP 13 EV

BTE and Thin Tube fitting solutions


The Thin Tube earmold is a great solution for individuals who are concerned with how an earmold appears cosmetically. The standard solution is appropriate for most hearing losses and ideal for severe to profound losses where acoustic seal is vital to the fitting's success.

Thin Tube Earmolds / Digital SLS / HS Silicone 60 / HS Silicone 40

FIND ME: Shop > Earmolds > Thin Tube Earmold

Tubing Lengths:

- 1
- 2
- 3, 3+
- 4, 4+
- 5, 5+

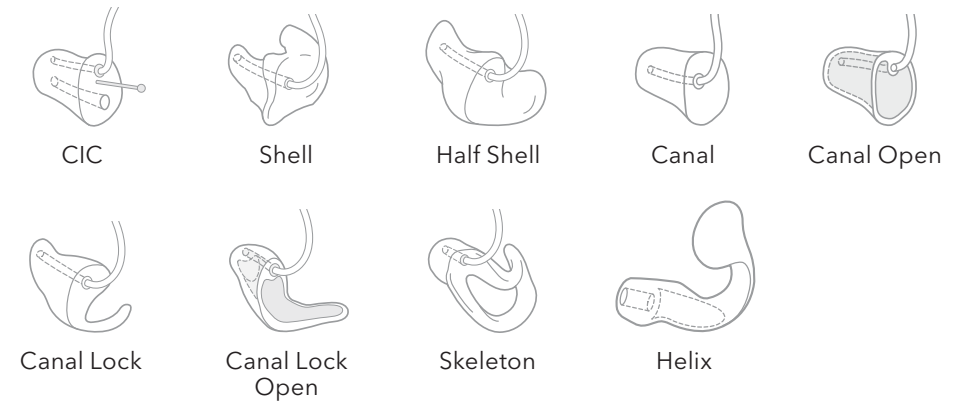


BTE Earmolds / Digital SLS / HS Silicone 60 / HS Silicone 40


FIND ME: Shop > Earmolds > Standard Earmold

Tubing:

- 13 Regular
- 13 Heavy
- 13 Extra Heavy
- 13 Dry Tubing
- 13 Regular Snap Connect
- 13 Heavy Snap Connect

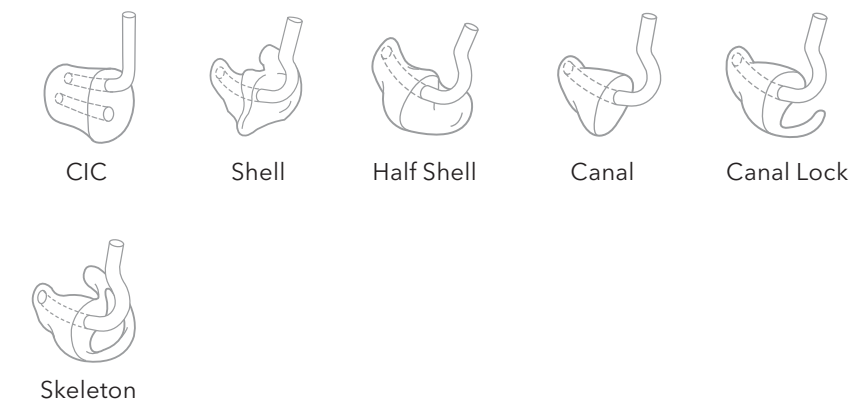

Stock solutions

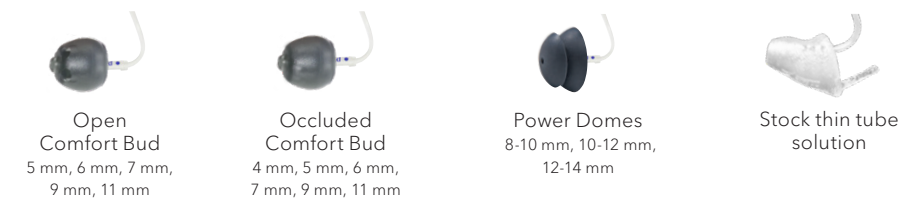

Open Comfort Bud
5 mm, 6 mm, 7 mm,
9 mm, 11 mm

Occluded Comfort Bud
4 mm, 5 mm, 6 mm,
7 mm, 9 mm, 11 mm

Power Domes
8-10 mm, 10-12 mm,
12-14 mm

Stock thin tube solution

Hearing aids do not restore natural hearing. Individual experiences vary depending on severity of hearing loss, accuracy of evaluation, proper fit and ability to adapt to amplification.

Information in this brochure is subject to change without notice.
©2022 Miracle-Ear, Inc. All rights reserved.

Miracle-Ear, Inc.
150 South Fifth Street
Suite 2300
Minneapolis, MN 55402



www.miracle-ear.com